

Thermo Scientific
SuperFrost and SuperFrost Printink

Color-coded archiving system

- Permanent lettering • Durable labeling
- Safe archiving of your slides and specimens

Permanent lettering and archiving of your slides

Based on the concept of the frosted writing area, Thermo Scientific™ SuperFrost® and Thermo Scientific™ SuperFrost Printink® slides have a colored writing area on one side. A wide range of colors (7 standard colors) allows several examination specimens to be labeled and archived simply and clearly.

- Time-consuming labeling with paper labels is no longer required
- Straightforward slide labeling
- The writing area is resistant to all common solvents
- Inkjet-printable writing area for SuperFrost Printink® slides
- Exact printing with one-dimensional or two-dimensional barcodes
- Ground slides are safe for use in automated laboratory systems
- Simplified workflow due to the various colors
- SuperFrost® slides may be sorted by color e.g. by patient, urgency, stain or examination method
- Simple recognition of specified characteristics such as year of examination or department, use etc.
- Adhesion or scratching of slides lying on top of each other is prevented

Product properties

- Manufactured from extra-white soda-lime glass with very low iron content
- Clean, flat and ready to use immediately
- Manufactured in accordance with ISO 8037/1 dimensions: 26 x 76 x 1.0 mm (+0/-1.0 mm)
- Tolerance for thickness: ± 0.05 mm
- Planarity under 30 μ m measured over 60 mm length



For the production of our slides, the highest-quality soda-lime raw materials are used to give first-rate, extra-white glass quality. A natural process of degradation may occur even when stored in accordance with the specifications. We therefore grant a warranty period of 12 months. The expiration date of the warranty period is printed on every box.




thermoscientific.com

© 2015 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries.

Slides & Specialty Glass

Thermo Fisher Scientific
Gerhard Menzel B.V. & Co. KG
 Saarbrückener Straße 248
 38116 Braunschweig, Germany

Tel. +49 (0) 531 59 00 80
 Fax +49 (0) 531 50 97 99
menzel.marketing@thermofisher.com

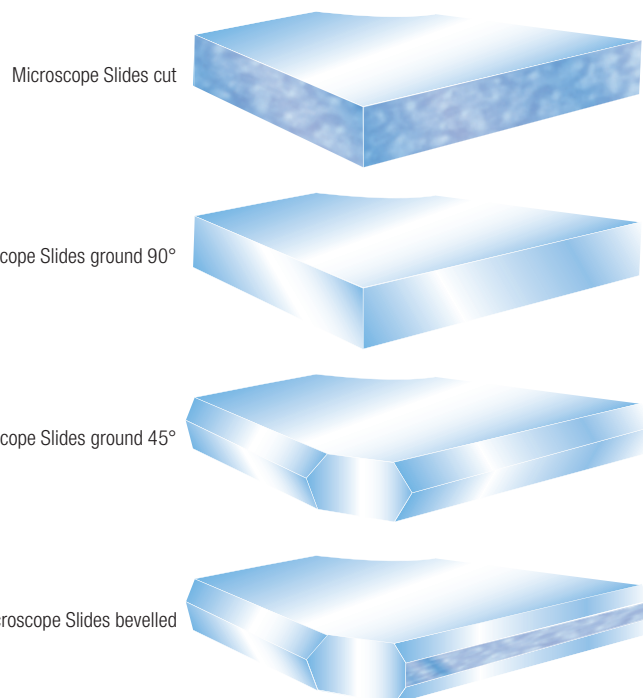

Stock list numbers:

Model

| | |
|---------------|--------------------|
| Cut | AAAA00008*##32E... |
| Ground 45° | ABAA00008*##32E... |
| Ground 90° | AGAA00008*##32E... |
| Ground 90° CC | AFAA00008*##32E... |
| Bevelled | ADAA00008*##32E... |

* 0 white 1 blue 2 pink 3 yellow 4 green 5 orange 6 violet

In order to protect the SuperFrost® slides from climatic influences, packaging suitable for tropical conditions is available for use in areas with high humidity. If this packaging is used, the delivered unit is vacuum packed. The glass is therefore protected from humidity and premature oxidation. This packaging extends the shelf-life up to 15 months. For further protection from external influences, the additional packaging variations "glass in cellophane" and/or "box in cellophane" are available.

Thermo
 S C I E N T I F I C

A Thermo Fisher Scientific Brand